


C₁-C₆ Alkanes and Alkenes

Column: GS-GasPro
15 m x 0.32 mm
J&W P/N: 113-4312
Carrier: Nitrogen at 25 cm/sec
Oven: 40°C for 1 min
40-200°C at 10°/min
Injector: Split 1:40, 100 µL, 200°C
Detector: FID, 200°C



1. Methane
2. Ethane
3. Ethylene
4. Propane
5. Propylene
6. *n*-Butane
7. 1-Butene
8. *n*-Pentane
9. 1-Pentene
10. *n*-Hexane
11. 1-Hexene

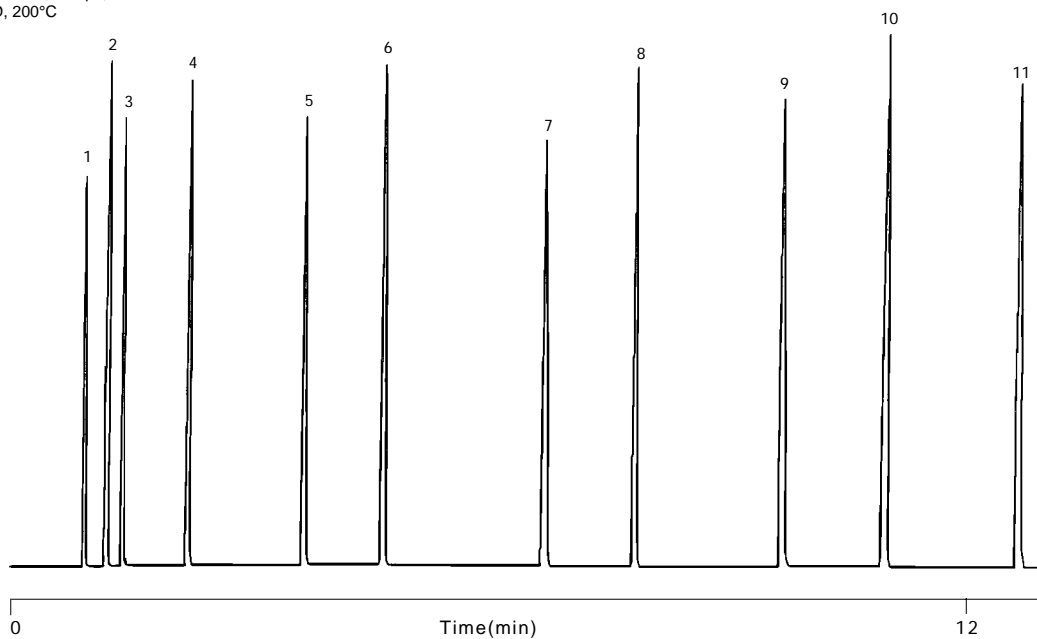

C533