


Aromatic Hydrocarbons

Column: **GS-GasPro**
30 m x 0.32 mm I.D.

J&W P/N: 113-4332
Carrier: Helium, 39.7 cm/s at 110°C
Oven: 110°C for 2 min
110-205°C at 10°/min
205-260°C at 25°C/min
260°C for 10 min
Injector: Split 1:20, 250°C
Detector: FID, 250°C

1. Pentane
2. Methylcyclopentane
3. 3-Methylpentane
4. Cyclohexane
5. Hexane
6. Heptane
7. Benzene
8. Isooctane
9. Octane
10. Toluene
11. Chlorobenzene
12. Nonane
13. Ethylbenzene
14. Styrene
15. *m,p*-Xylene
16. Decane
17. *o*-Xylene
18. Cumene
19. Propylbenzene
20. 1,2,4-Trimethylbenzene
21. 1,3,5-Trimethylbenzene
22. Undecane

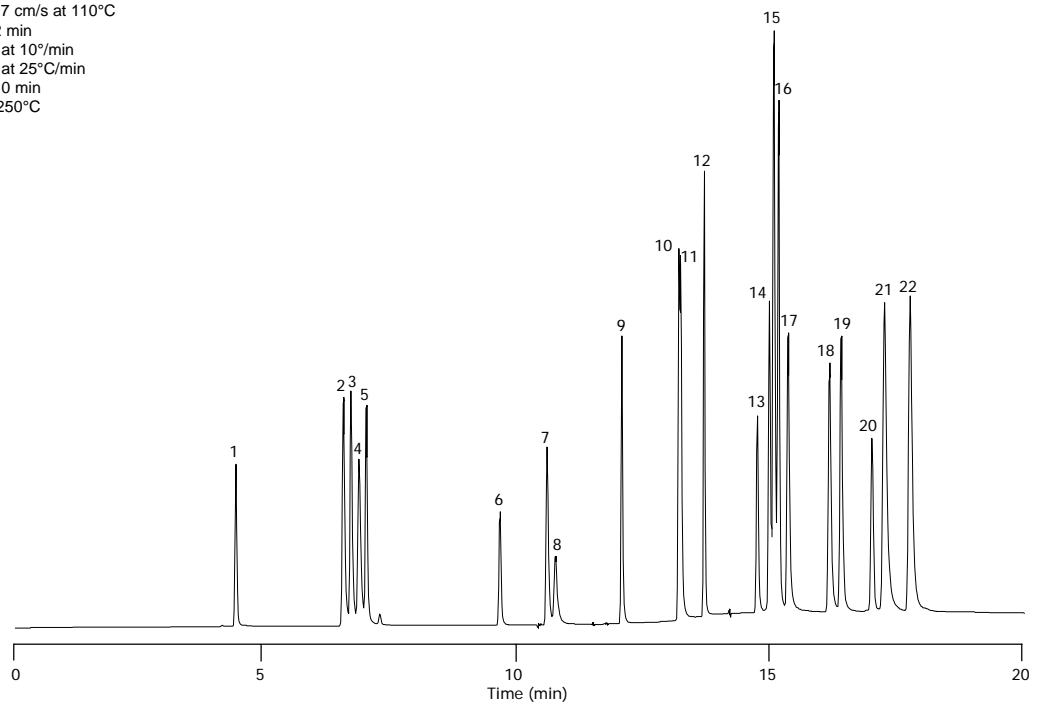

C676