


Ethanolamines Rtx®-5 Amine

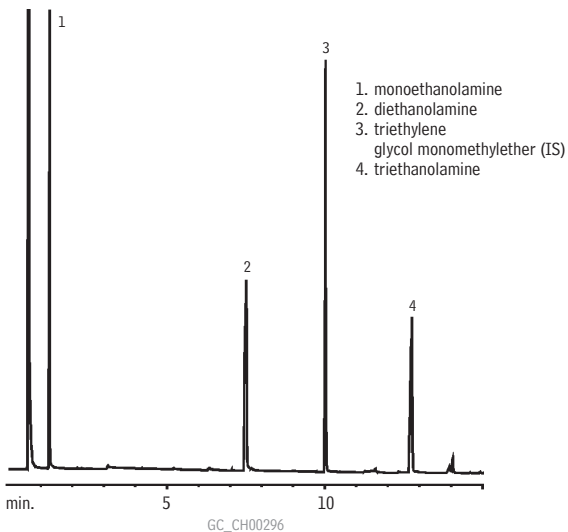

15m, 0.25mm ID, 0.50 μ m Rtx®-5 Amine (cat.# 12335)
1.0 μ L split injection of ethanolamine mix in methanol
On-column conc.: 34ng

Oven temp.: 50°C (hold 2 min.) to 180°C @
10°C/min. (hold 2 min.)
Inj. & det. temp.: 280°C / 300°C
Carrier gas: hydrogen
Linear velocity: 43cm/sec. set @ 50°C
FID sensitivity: 6.4 x 10⁻¹¹ AFS
Split ratio: 58:1